Aurora B IHC Antibody

Autora Bille Autobady					
Rabbit Polyclonal Antigen Affinity Purified Catalog No. IHC-00003-T		Protein IE T GeneID	D Q96GD4 9212		
Lot No. IHC-00003-T-1		T-1		BEIHYL	
				LABOHATOHTES, INC	
APPLICATIONS		, ICC-IF			
SPECIES REACTI		Human			
AMOUNT					
CONCENTRATIC					
STORAGE/SHELI		2 – 8° C / 1 year from date of receipt			
PHYSICAL STAT	E Liqu	Liquid			
BUFFER	Tris	Tris-buffered Saline containing 0.1% BSA and 0.09% Sodium Azide			
ISOTYPE	lgG	IgG			
ORIGIN	USA	USA			
PRODUCTION	Anti	Antibody was affinity purified using an epitope specific to Aurora B immobilized on solid support.			
PROCEDURES	The epitope recognized by IHC-00003-T maps to a region between residue 300 and the terminus (residue 344) of human Aurora B using the numbering given in SwissProt entry Q96GD4 (GeneID 9212).				
APPLICATIONS		Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.			
	Imr	nunohistochemistry	1:100 - 1:500		
	lmr (ICC		1:100 - 1:500		
APPLICATION N	Epit	Formalin-fixed Paraffin-Embedded Tissue 1:100 – 1:500: Epitope exposure is recommended. Epitope exposure with citrate buffer will enhance staining.			
		nalin-fixed cells 1:100 neabilization is recomn			
	Like	ly to work with frozen	sections.		
	In se	ome cases, the antibod	y may be diluted furth	er than indicated.	
ADDITIONAL IN		s://www.bethyl.com/p the link above to view		itations, and other product specific information.	


This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc. Eric McIntush, PhD | Chief Scientific Officer Date: June 21, 2019

Bethyl Laboratories, Inc. • 25043 West FM 1097 • Montgomery, TX 77356 • 800.338.9579 • 936.597.6111 • 866.597.6105 (FAX) • www.bethyl.com • technical@bethyl.com

Aurora B IHC Antibody



Detection of human Aurora B by immunohistochemistry. *Sample:* FFPE section of human colon adenocarcinoma. *Antibody:* Affinity purified rabbit anti-Aurora B (Cat. No. IHC-00003-T) used at a dilution of 1:100. *Detection:* DAB

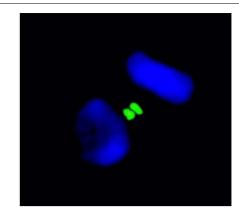

Detection of human Aurora B by immunocytochemistry. *Sample:* NBF-fixed HeLa Cells. *Antibody:* Affinity purified rabbit anti-Aurora B (Cat. No. IHC-00003-T) used at a dilution of 1:100. *Detection:* goat anti-rabbit IgG (h&l) FITC (Cat. No. A120-101F) used at a dilution of 1:100.