

# Filamin B IHC Antibody

Rabbit Polyclonal

Antigen Affinity Purified

Protein ID NP\_001448.2

Catalog No. IHC-00355-T

GeneID 2317

Lot No. IHC-00355-T-1



**APPLICATIONS** IHC

**SPECIES REACTIVITY** Human

**AMOUNT** 10 µl

**CONCENTRATION** 100 µg/ml

**STORAGE/SHELF LIFE** 2 – 8° C / 1 year from date of receipt

**PHYSICAL STATE** Liquid

**BUFFER** Tris-buffered Saline containing 0.1% BSA and 0.09% Sodium Azide

**ISOTYPE** IgG

**ORIGIN** USA

**PRODUCTION PROCEDURES** Antibody was affinity purified using an epitope specific to Filamin B immobilized on solid support.

The epitope recognized by IHC-00355-T maps to a region between residue 2552 and 2602 of human filamin B, beta (actin binding protein 278) using the numbering given in entry NP\_001448.2 (GeneID 2317).

Antibody concentration was determined by extinction coefficient: absorbance at 280 nm of 1.4 equals 1.0 mg of IgG.

**APPLICATIONS** Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.

Immunohistochemistry 1:100 – 1:500

**APPLICATION NOTES** Epitope exposure is recommended.

Epitope exposure with citrate buffer will enhance staining.

Likely to work with frozen sections.


In some cases, the antibody may be diluted further than indicated.

**IHC HUMAN CONTROLS** Colon Carcinoma, Prostate Carcinoma


**ADDITIONAL INFO** <https://www.bethyl.com/product/IHC-00355-T>

Use the link above to view SDS, a current list of citations, and other product specific information.

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc.  
Eric McIntush, PhD | Chief Scientific Officer Date: June 21, 2019

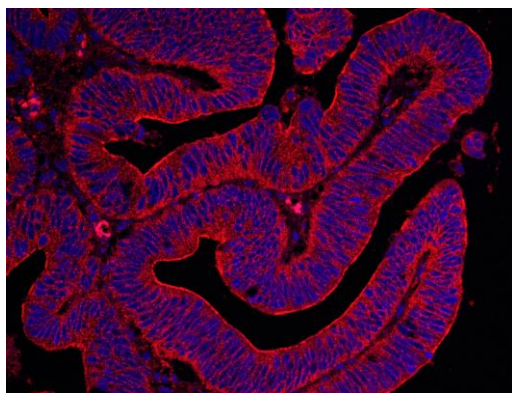

**Detection of human Filamin B by immunohistochemistry.**  
*Sample:* FFPE section of human prostate carcinoma.  
*Antibody:* Affinity purified rabbit anti-Filamin B (Cat. No. IHC-00355-T) used at a dilution of 1:100. *Detection:* Red-fluorescent goat anti-rabbit IgG highly cross-adsorbed Antibody Hilyte Plus™ 555 used at a dilution of 1:100.

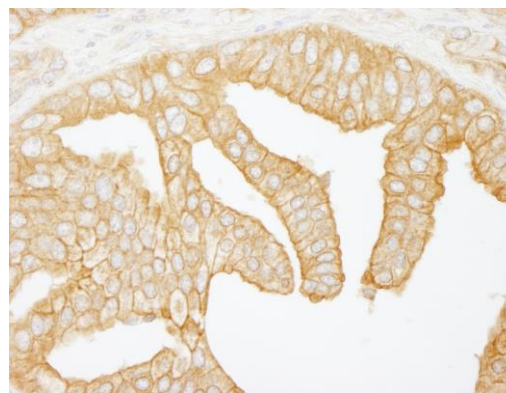

**Detection of human Filamin B by immunohistochemistry.**  
*Sample:* FFPE section of human prostate carcinoma.  
*Antibody:* Affinity purified rabbit anti-Filamin B (Cat. No. IHC-00355-T) used at a dilution of 1:500. *Detection:* DAB